


Bacterial Cell Membrane FAMES on Rxi®-17Sil MS

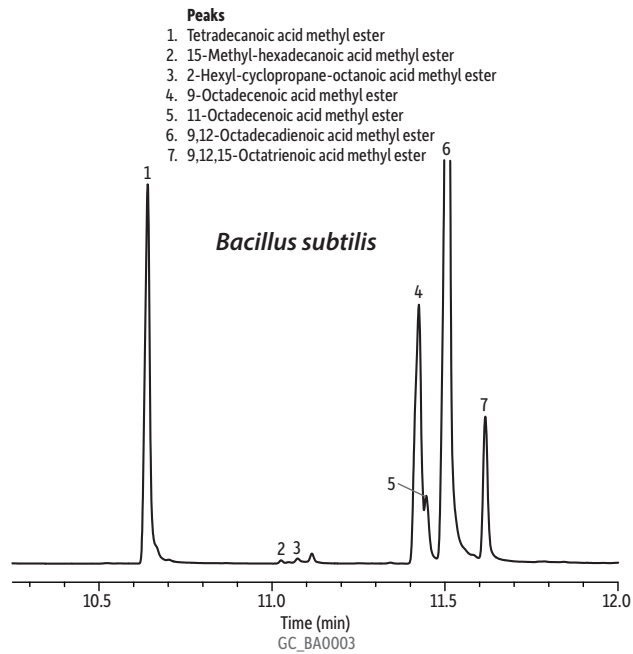

Column Rxi®-17Sil MS, 20 m, 0.15 mm ID, 0.15 µm (cat.# 43821)
Sample FAMES in bacterial cell wall lecithin
Injection
Inj. Vol.: 1 µL split (split ratio 100:1)
Inj. Temp.: 250 °C
Oven
Oven Temp.: 70 °C (hold 4 min) to 250 °C at 25 °C/min (hold 10 min)
Carrier Gas He, constant linear velocity
Linear Velocity: 50 cm/sec
Detector DANI Master GC-TOFMS @ 25 Hz
Acknowledgement Data provided by Dr. Radomír Čabala, Head of the Toxicology Department at the General University Hospital in Prague, Czech Republic

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
1-814-353-1300 • 1-800-356-1688 • Fax: 1-814-353-1309 • www.restek.com