


Column Test Mix  
Rxi®-5Sil MS

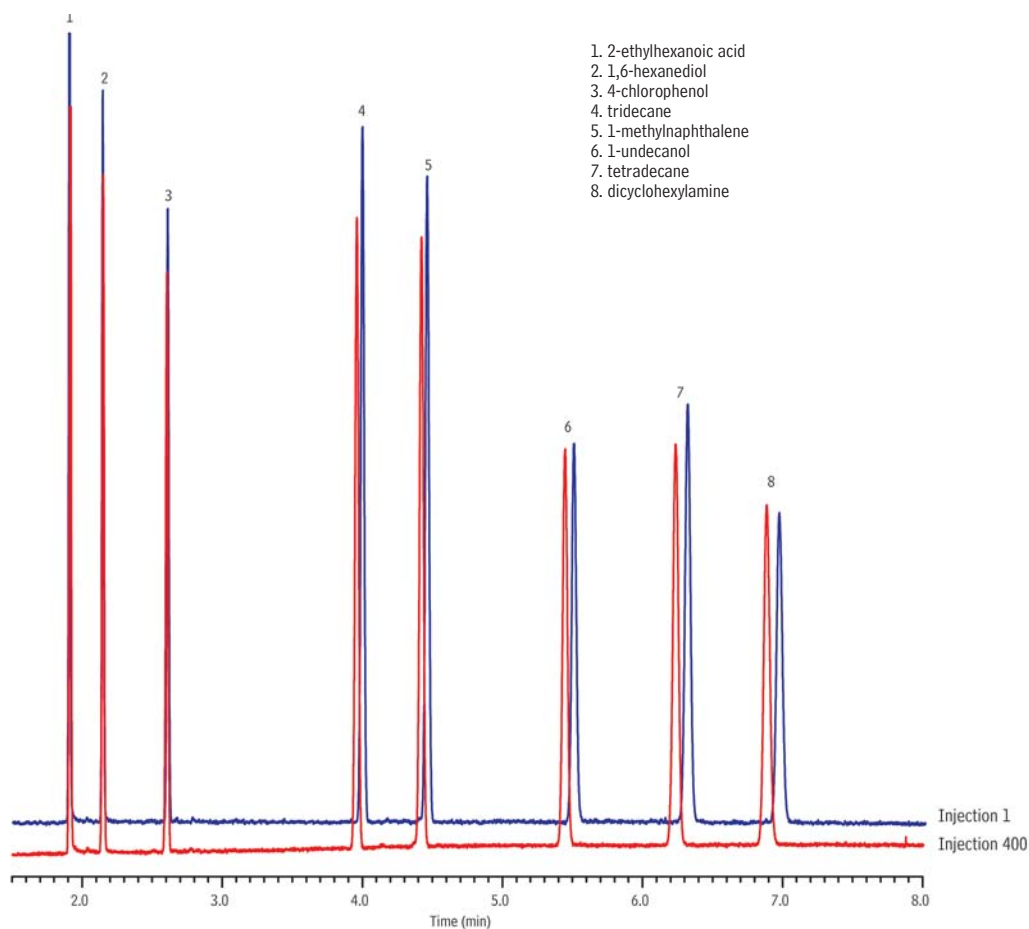

GC\_CF01131

Column: Rxi®-5Sil MS, 30m, 0.25mm ID, 0.25 $\mu$ m (cat.# 13623)  
Sample: Column Test Mix (cat.# 35226)  
Inj.: 1.0 $\mu$ L split (split ratio 1:60), 4mm recessed gooseneck liner (cat.# 20983)  
Inj. temp.: 250°C  
Carrier gas: helium, constant pressure  
Linear velocity: 36cm/sec @ 125°C  
Oven temp.: 125°C  
Det: FID @ 320°C  
Instrument: Agilent 6890