


Alcohols Rtx®-5

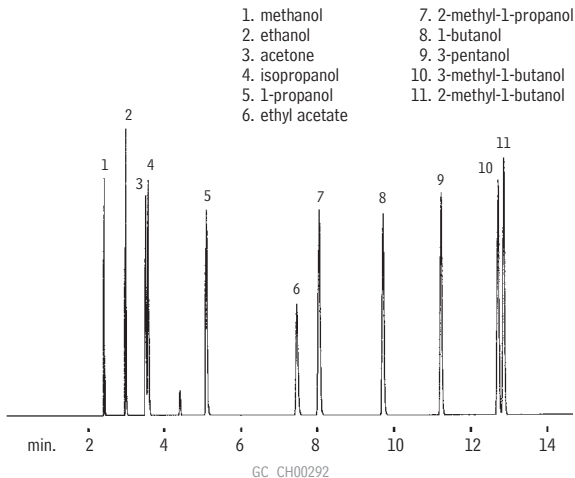

60m, 0.32mm ID, 1.0 μ m Rtx®-5 (cat.# 10257)
0.03 μ L split injection of a neat alcohols mix

Oven temp.: 25°C (hold 4 min.) to 80°C @
8°C/min. (hold 5 min.)
Inj. & det. temp.: 200°C
Carrier gas: hydrogen
Linear velocity: 40cm/sec. (flow rate: 3.2cc/min.)
FID sensitivity: 128 x 10¹¹ AFS
Split ratio: 40:1