


Hexamethylenediamine Stabilwax®-DB

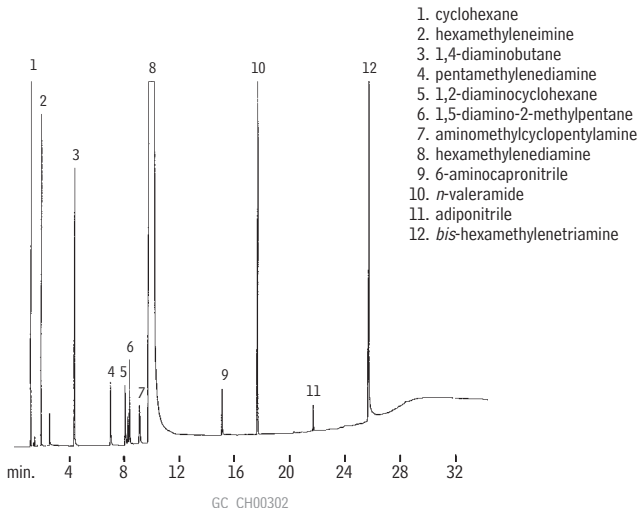

30m, 0.32mm ID, 0.25 μ m Stabilwax-DB (cat.# 10824)
0.4 μ L direct injection of a neat hexamethylenediamine sample
On-column conc.: 10 to 1000ng/component

Oven temp.: 95°C (hold 6 min.) to 235°C @
7°C/min. (hold 4 min.)

Inj. & det. temp.: 250°C
Carrier gas: hydrogen
Linear velocity: 40cm/sec.
FID sensitivity: 2 x 10⁻¹¹ AFS