


# Nitrosamines Stabilwax®-DB

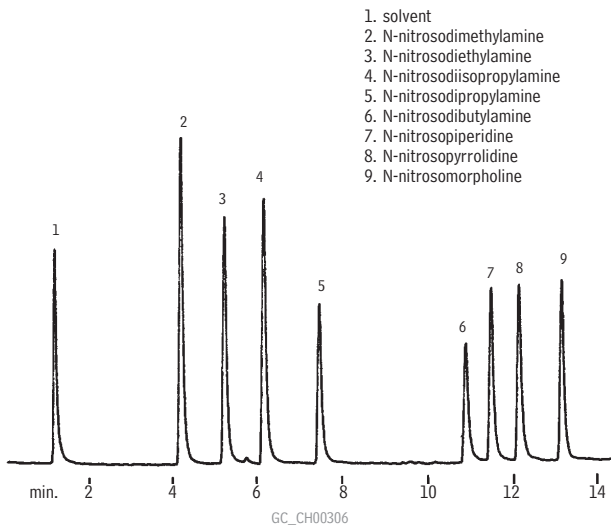

1. solvent
2. N-nitrosodimethylamine
3. N-nitrosodiethylamine
4. N-nitrosodiisopropylamine
5. N-nitrosodipropylamine
6. N-nitrosodibutylamine
7. N-nitrosopiperidine
8. N-nitrosopyrrolidine
9. N-nitrosomorpholine

60m, 0.53mm ID, 1.0 $\mu$ m Stabilwax®-DB (cat.# 10858)  
Direct injection of nitrosamines  
Conc.: 1.0 $\mu$ g/mL

Oven temp.: 100°C (hold 1 min.) to 170°C @ 5°C/min.  
Inj. & det. temp.: 200°C  
Carrier gas: helium  
Linear velocity: 100cm/sec. (flow rate: 15cc/min.)  
Detector: TSD