


## Amines/Alcohols/Chlorides Stabilwax®

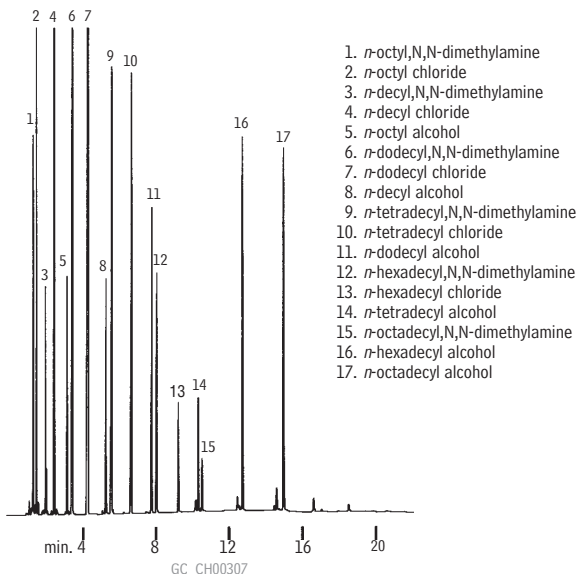

30m, 0.53mm ID, .25 $\mu$ m Stabilwax® (cat.# 10625)  
0.5 $\mu$ L split injection

Oven temp.: 100°C to 250°C @  
8°C/min. (hold 5 min.)  
Inj. & det. temp.: 250°C  
Carrier gas: hydrogen  
Linear velocity: 40cm/sec. (flow rate: 5.2cc/min.)  
FID sensitivity: 128 x 10<sup>11</sup> AFS  
Split ratio: 40:1

Restek Corporation 110 Benner Circle Bellefonte, PA 16823  
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com