


Chlorofluorocarbons Rtx[®]-1

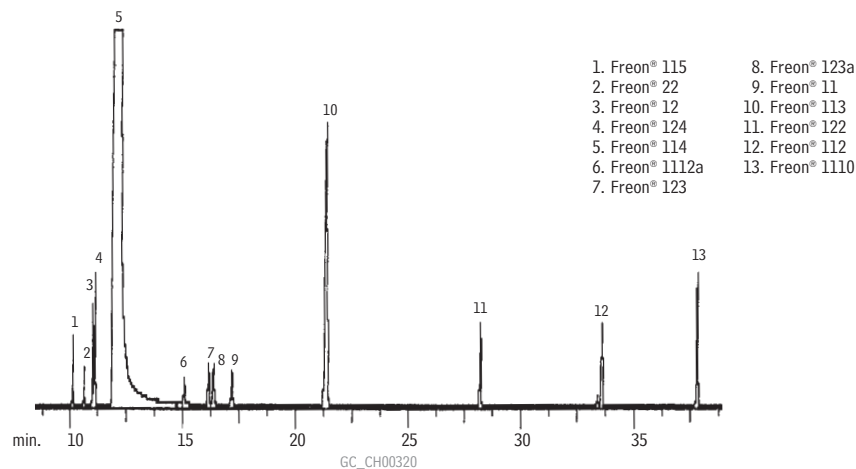

105m, 0.25mm ID, 1.0 μ m Rtx[®]-1 (cat.# 10159)
1.0 μ L split injection of a 99.9% pure Freon[®] 114 standard.

Oven temp.: 35°C (hold 20 min.) to 220°C
@ 20°C/min. (hold 30 min.)
Inj. & det. temp.: 250°C
Det.: full scan
Carrier gas: helium
Linear velocity: 20cm/sec. set @ 35°C
Split vent: 30mL/min.