


# Freon® & Ethylene Oxide Rtx®-1

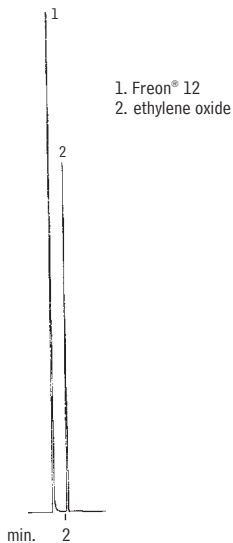

GC\_CH00324

30m, 0.53mm ID, 5.0 $\mu$ m Rtx®-1 (cat.# 10179)  
50 $\mu$ L split injection of Freon® 12 and ETO

Oven temp.:	25°C
Inj. & det. temp.:	290°C
Carrier gas:	hydrogen
Linear velocity:	40cm/sec.
FID sensitivity:	4 x 10 <sup>-11</sup> AFS
Split vent:	40:1

Restek Corporation 110 Benner Circle Bellefonte, PA 16823  
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com