


Solvents HayeSep® Q

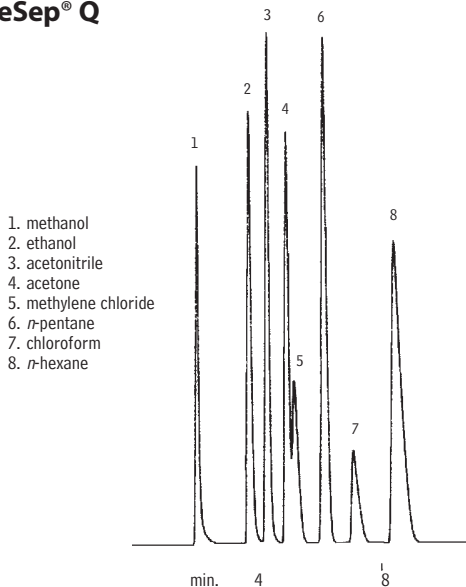

1. methanol
2. ethanol
3. acetonitrile
4. acetone
5. methylene chloride
6. *n*-pentane
7. chloroform
8. *n*-hexane

GC_CH00328

2m, 1mm ID HayeSep® Q (cat.# 19017)
1µL direct injection of a neat solvent mixture

Oven temp.: 80°C to 180°C @
16°C/min. (hold 5 min.)
Inj. & det. temp.: 200°C
Carrier gas: helium
Flow: 20mL/min. set @ 40°C
FID sensitivity: 512 x 10⁻¹¹ AFS

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com