



Pentane (lot purity analysis) MXT®-1

LOT A

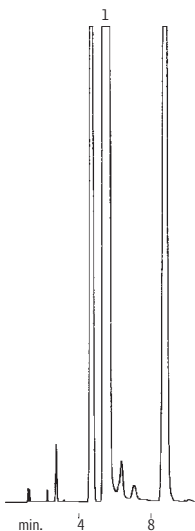

LOT B

1. pentane

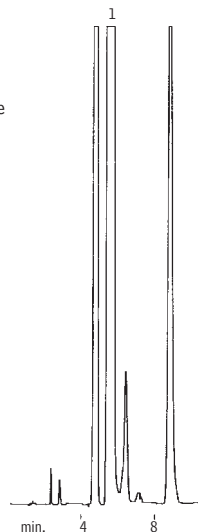

GC_CH00332

60m, 0.53mm ID, 5.0 μ m MXT®-1 (cat.# 70183)
1.0 μ L direct injection of high purity pentane

Oven temp.: 40°C
Inj. & det. temp.: 200°C
Carrier gas: helium
Linear velocity: 40cm/sec. (flow rate: 10cc/min.)
FID sensitivity: 4 x 10¹¹ AFS

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com