


Substituted Pyridine Stabilwax®

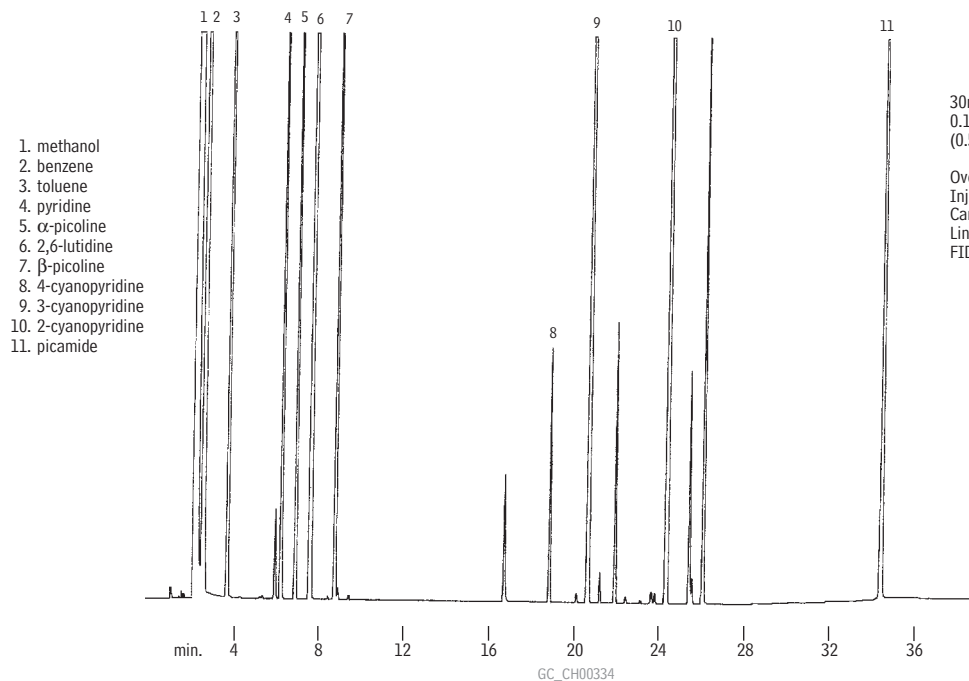

30m, 0.53mm ID, 1.0 μ m Stabilwax® (cat.# 10655)
0.1 μ L direct injection of a substituted pyridine standard
(0.5-5% per component), Uniliner® liner

Oven temp.: 60°C to 225°C @ 5°C/min.
Inj. & det. temp.: 230°C
Carrier gas: hydrogen
Linear velocity: 40cm/sec. (flow rate: 5.1cc/min.)
FID sensitivity: 32 x 10⁻¹¹ AFS