


Volatile Solvents Rtx[®]-Wax

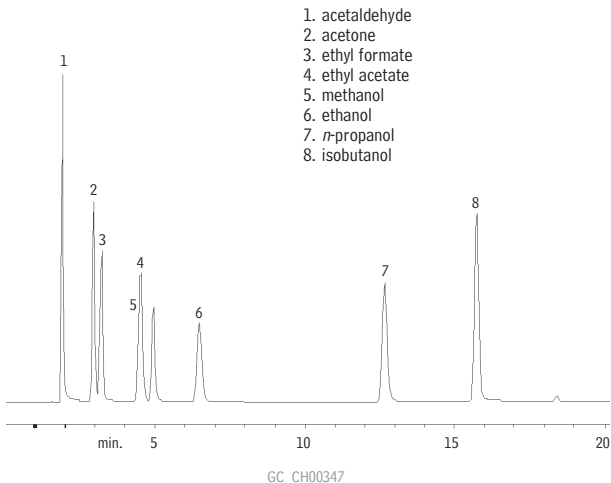

30m, 0.53mm ID, 1.0 μ m Rtx[®]-Wax (cat.# 12455)
0.1 μ L split injection of volatile solvents

Oven temp.: 20°C (hold 10 min.) to 100°C @ 5°C/min.
Inj. & det. temp.: 200°C
Carrier gas: helium
Linear velocity: 40cm/sec. set @ 20°C
FID sensitivity: 8 x 10⁻¹¹ AFS
Split ratio: 10:1