


Solvents, Industrial Stabilwax®-DA

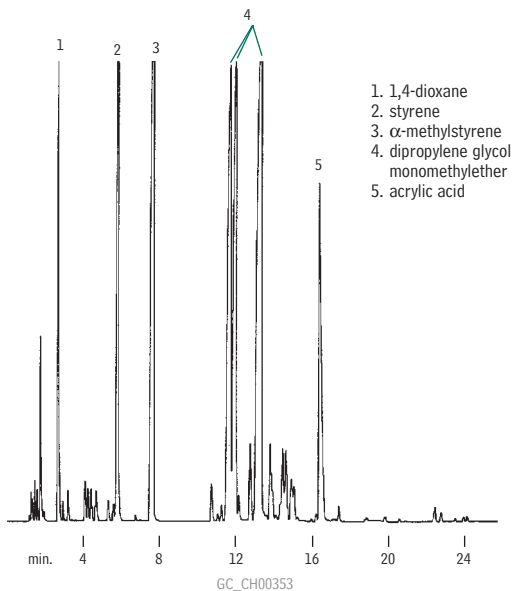

1. 1,4-dioxane
2. styrene
3. α -methylstyrene
4. dipropylene glycol monomethylether
5. acrylic acid

30m, 0.53mm ID, 1.0 μ m Stabilwax®-DA (cat.# 11055)
0.1 μ L direct injection of industrial solvents

Oven temp.: 95°C (hold 4 min.) to 250°C @ 4°C/min.
Inj. & det. temp.: 250°C
Carrier gas: helium
Linear velocity: 40cm/sec. (flow rate: 5cc/min.)
FID sensitivity: 32 x 10⁻¹¹ AFS