


Styrene Impurities Stabilwax®

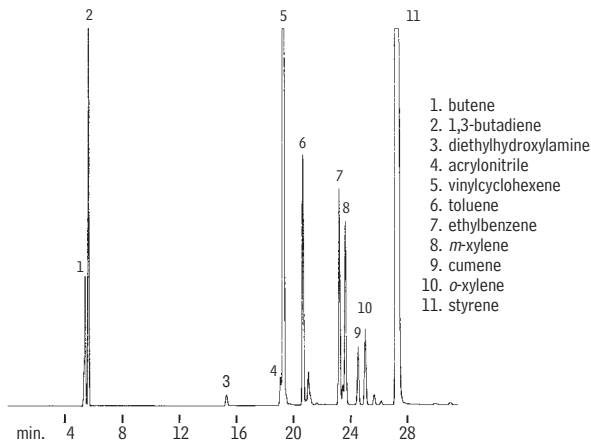

GC_CH00356

30m, 0.53mm ID, 0.50 μ m Stabilwax® (cat.# 10640)
0.5 μ L split injection of a 95% pure styrene

Oven temp.: 40°C (hold 10 min.) to 150°C
@ 10°C/min. (hold 15 min.)
Inj. & det. temp.: 150°C
Carrier gas: helium
Linear velocity: 20cm/sec. set @ 40°C
FID sensitivity: 16 x 10⁻¹¹ AFS
Split vent: 40cc/min.

Permission to publish this chromatogram granted by
Copolymer Rubber and Chemical Corp.

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com