


Silanes Rtx®-200

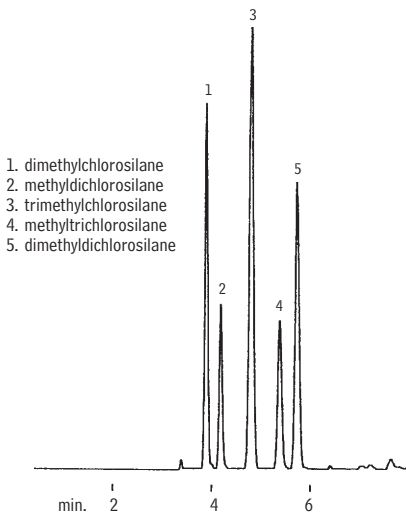

GC_CH00362

60m, 0.53mm ID, 3.0 μ m Rtx®-200 (cat.# 15088)
0.5 μ L split injection of methyl silanes

Oven temp.: 40°C to 250°C @ 8°C/min.
(hold 5 min.)
Inj. / det. temp.: 200°C / 270°C
Carrier gas: hydrogen
Linear velocity: 40cm/sec.
FID sensitivity: 1.02 x 10⁹ AFS
Split vent: 40cc/min.

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com