


**Uber Test Mix
Rxi®-5Sil MS**

Column: Rxi®-5Sil MS, 30m, 0.25mm ID, 0.25 μ m (cat.# 13623)
Sample: Uber Test Mix
Inj.: 1.0 μ L, split (split vent flow 200mL/min.), 4mm gooseneck splitless liner (cat.# 20798)
Inj. temp.: 250°C
Carrier gas: hydrogen
Linear velocity: 40cm/sec. @ 65°C
Oven temp.: 65°C isothermal
Det: FID @ 325°C

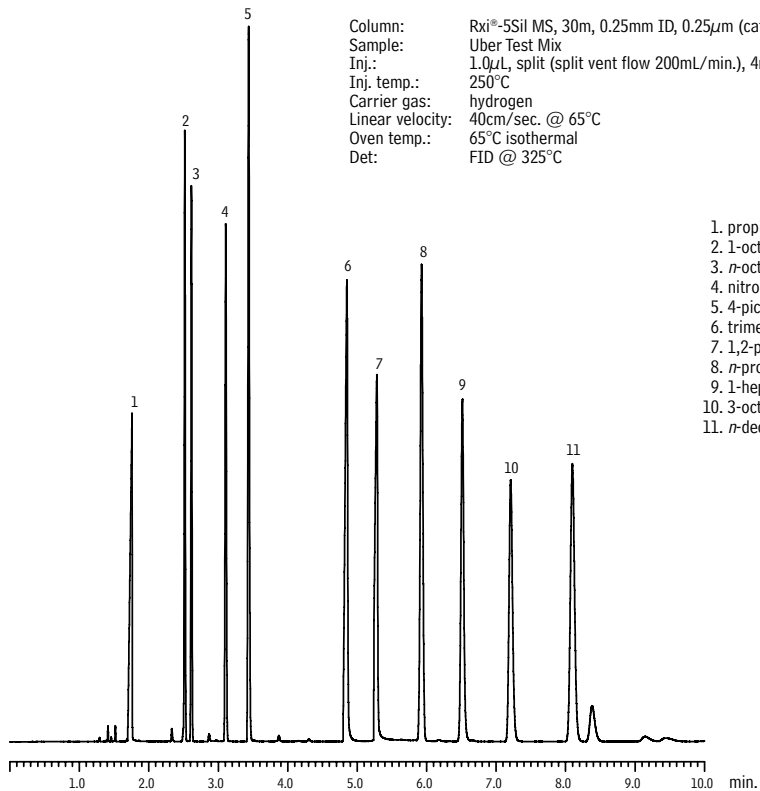

1. propionic acid
2. 1-octene
3. *n*-octane
4. nitrobutane
5. 4-picoline
6. trimethyl phosphate
7. 1,2-pentanediol
8. *n*-propylbenzene
9. 1-heptanol
10. 3-octanone
11. *n*-decane

GC_CH01097