


Nitroaromatics

US EPA Method 609

Rtx®-5

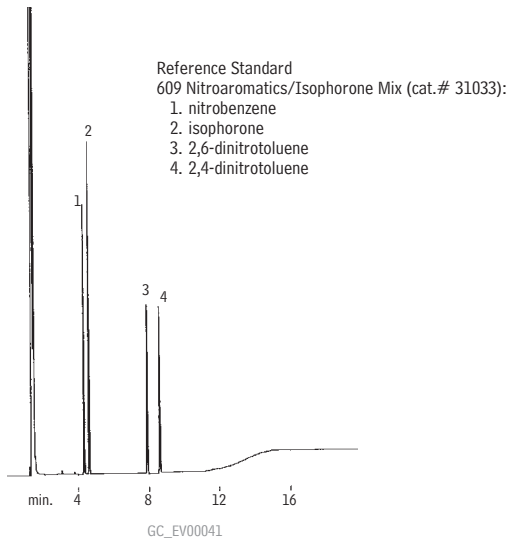

15m, 0.53mm ID, 1.50 μ m Rtx®-5 (cat.# 10267)

0.25 μ L split injection of nitroaromatics.

On-column concentration 28ng/component.

Oven temp.:	100°C to 300°C @ 15°C/min. (hold 10 min.)
Inj. / det. temp.:	250°C / 300°C
Carrier gas:	helium
Linear velocity:	23cm/sec. set @ 40°C
FID sensitivity:	range = 4, atten. = 3

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com