


DRO Mix Rtx®-5

30m, 0.53mm ID, 1.0 μ m Rtx®-5 (cat.# 10255)
Direct injection of DRO Mix (cat.# 31064)
Conc.: 50ng/ μ L

Oven temp.: 40°C (hold 2 min.) to 310°C
@ 12°C/min. (hold 10 min.)
Inj. / det. temp.: 280°C / 310°C
Carrier gas: helium @ 65cm/sec.
FID sensitivity: 8 x 10⁻¹¹ AFS

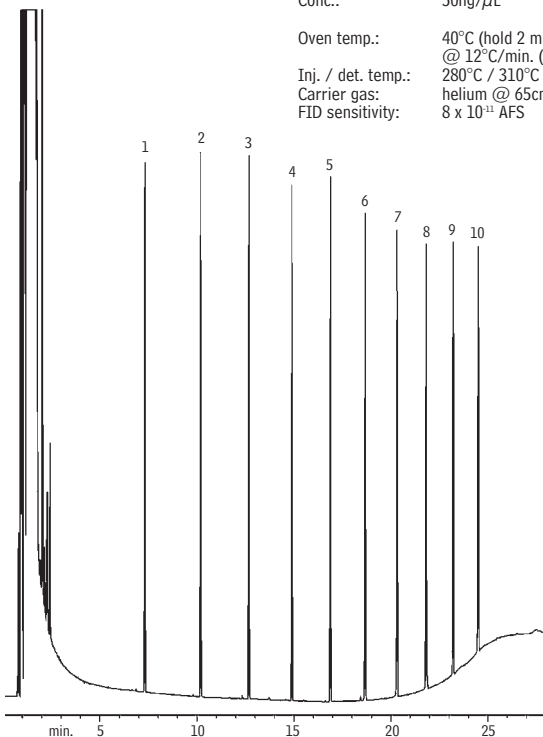

1. C10
2. C12
3. C14
4. C16
5. C18
6. C20
7. C22
8. C24
9. C26
10. C28

GC_EV00069