


Nitrogen-Containing Herbicides Rtx[®]-5

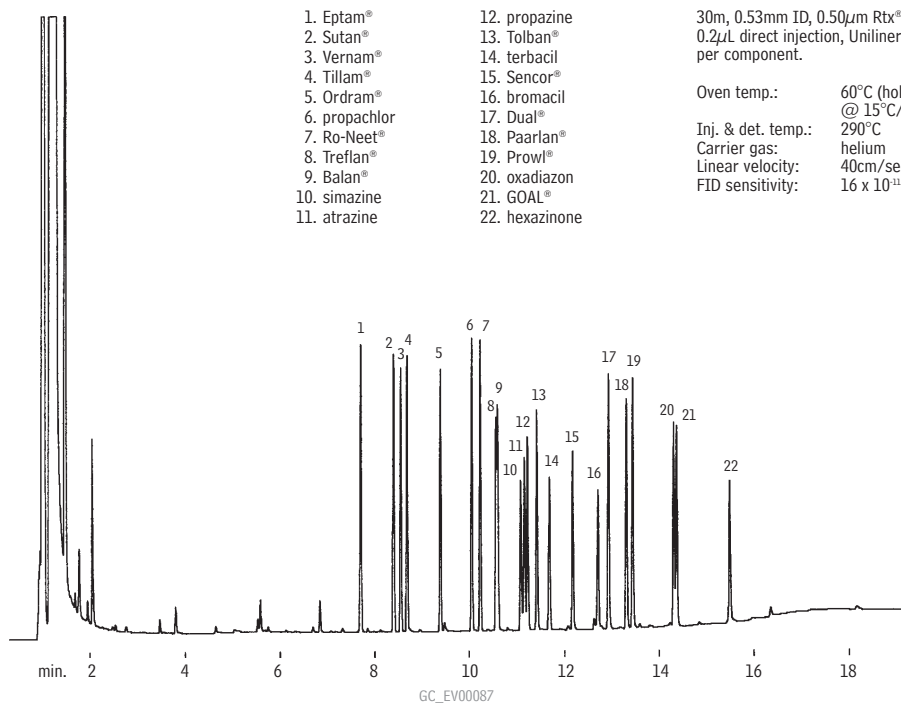

30m, 0.53mm ID, 0.50 μ m Rtx[®]-5 (cat.# 10240)
0.2 μ L direct injection, Uniliner[®] inlet liner, concentration approximately 10ng per component.

Oven temp.: 60°C (hold 1 min.) to 290°C
@ 15°C/min. (hold 5 min.)
Inj. & det. temp.: 290°C
Carrier gas: helium
Linear velocity: 40cm/sec. (flow rate: 5.2cc/min.)
FID sensitivity: 16 x 10¹¹ AFS