


Nitrogen-Containing Herbicides Rtx®-35

30m, 0.53mm ID, 0.50µm Rtx®-35 (cat.# 10440)
0.2µL direct injection, Uniliner® liner, concentration approximately 10ng per component.

Oven temp.: 60°C (hold 1 min.) to 290°C
@ 15°C/min. (hold 5 min.)
Inj. & det. temp.: 290°C
Carrier gas: helium
Linear velocity: 40cm/sec. (flow rate: 5.2cc/min.)
FID sensitivity: 16 x 10⁻¹¹ AFS

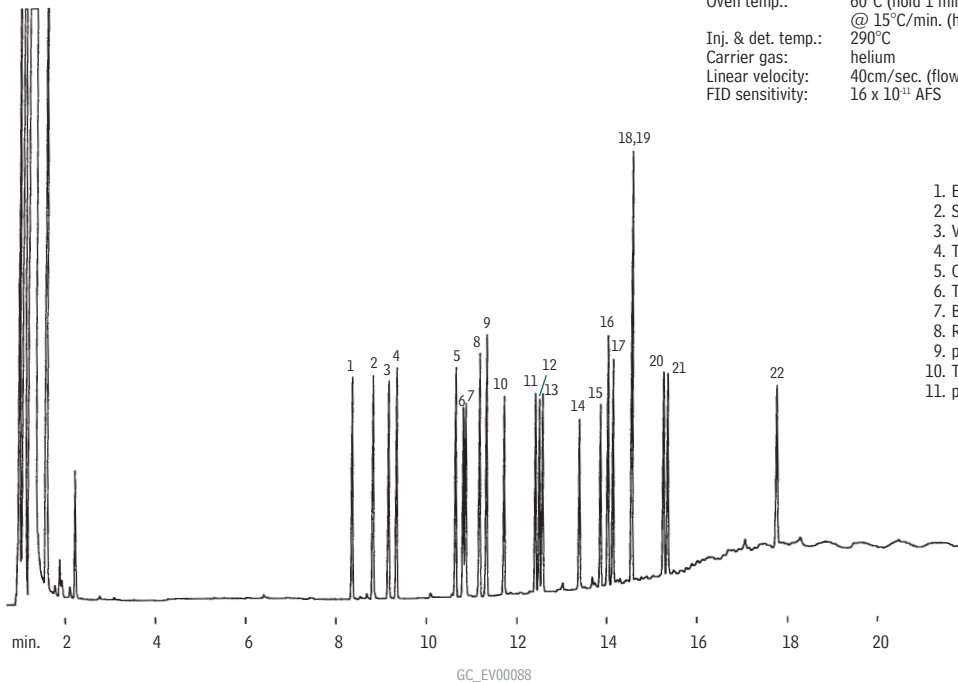

- | | |
|---------------|----------------|
| 1. Eptam® | 12. atrazine |
| 2. Sutan® | 13. simazine |
| 3. Vernam® | 14. terbacil |
| 4. Tillam® | 15. Sencor® |
| 5. Ordram® | 16. Dual® |
| 6. Treflan® | 17. Paarlan® |
| 7. Balan® | 18. Prowl® |
| 8. Ro-Neet® | 19. bromacil |
| 9. propachlor | 20. oxadiazol |
| 10. Tolban® | 21. GOAL® |
| 11. propazine | 22. hexazinone |