


Recovery of Chlorinated Pesticides from Solvent Extracts Using CarboPrep™ 90

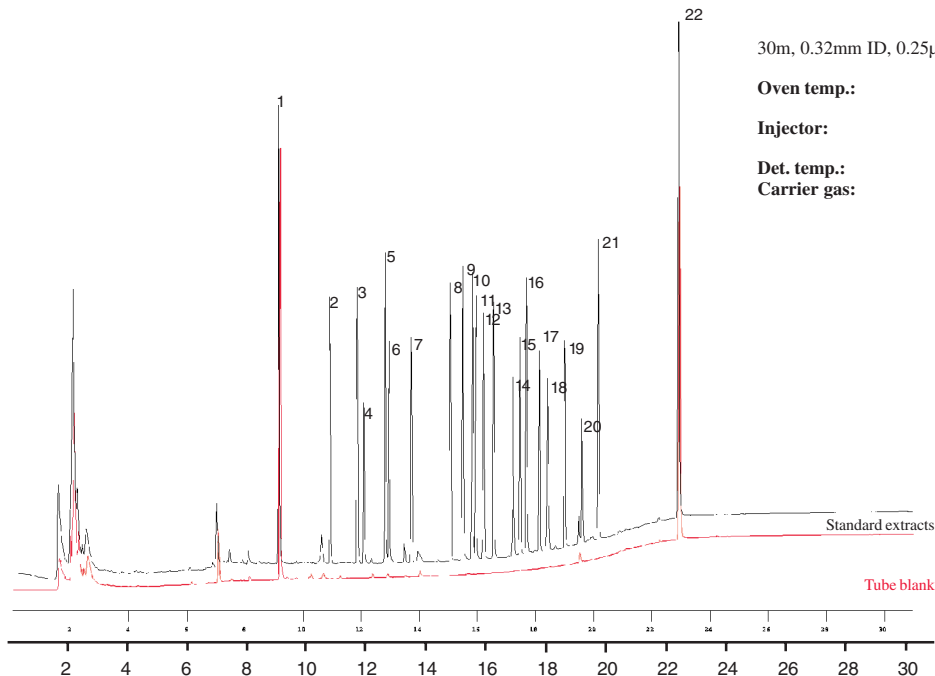

30m, 0.32mm ID, 0.25µm Rtx® CLPesticides2 (cat.# 11324).

Oven temp.: 120°C (hold 1 min.) to 300°C @
9°C/min. (hold 10 min.)

Injector: direct injection using Uniliner® sleeve
(cat.# 20335), 210°C.

Det. temp.: ECD 310°C

Carrier gas: helium

1. 2,4,5,6-tetrachloro-m-xylene
2. α-BHC
3. γ-BHC
4. β-BHC
5. δ-BHC
6. heptachlor
7. aldrin
8. heptachlor epoxide
9. γ-chlordane
10. α-chlordane
11. 4,4'-DDE
12. endosulfan I
13. dieldrin
14. endrin
15. 4,4'-DDD
16. endosulfan II
17. 4,4'-DDT
18. endrin aldehyde
19. methoxychlor
20. endosulfan sulfate
21. endrin ketone
22. decachlorobiphenyl