


# Lemon Oil

## Rt- $\beta$ DEXsm™

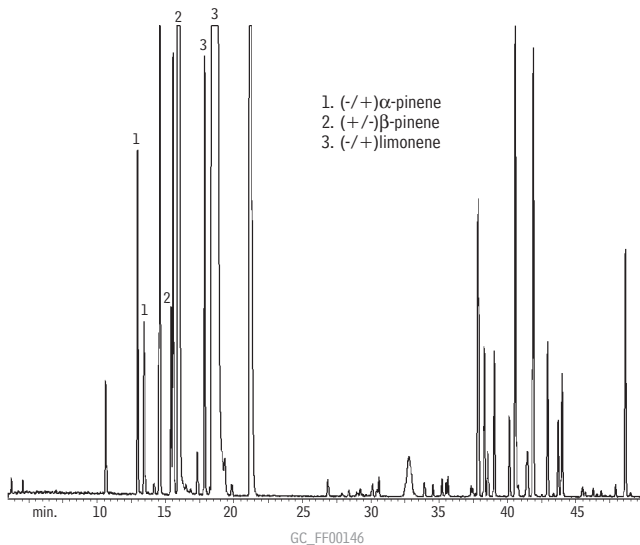

30m, 0.32mm ID, 0.25 $\mu$ m  
Rt- $\beta$ DEXsm™ (cat.# 13104)  
Oven temp.: 40°C (hold 1 min.) to  
200°C @ 2°C/min. (hold 3 min.)  
Carrier gas: hydrogen  
Linear velocity: 80cm/sec.  
Det.: FID set @ 220°C

Restek Corporation 110 Benner Circle Bellefonte, PA 16823  
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com