


# FAMEs (Black Currant Seed Oil)

## FAMEWAX™

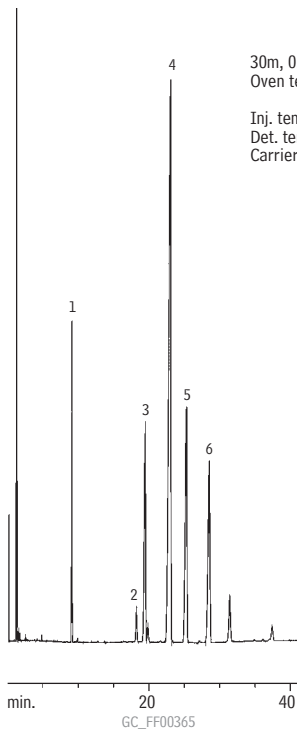

30m, 0.25mm ID, 0.25 $\mu$ m FAMEWAX™ (cat.# 12497).  
Oven temp.: 165°C (hold 30 min.) to 220°C @  
1.5°C/min. (hold 15 min.)  
Inj. temp.: 225°C  
Det. temp.: 230°C  
Carrier gas: helium @ 40cm/sec.

1. C16:0
2. C18:0
3. C18:1n9
4. C18:2n6
5. C18:3n6
6. C18:3n3