



Flavor & Fragrance Compounds Rt-CW20M™ F&F

Carbowax® column

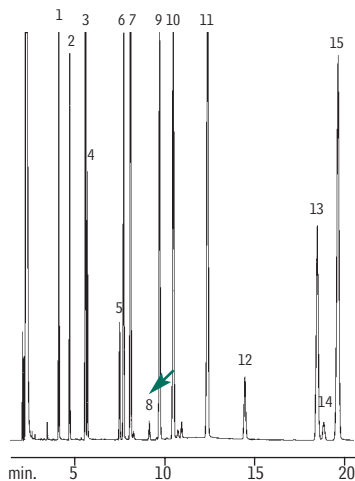

Rt-CW20M™ F&F column

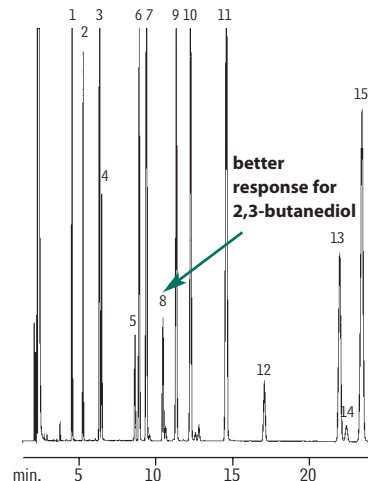

1. methyl heptanoate
2. hexanol
3. methyl octanoate
4. nonanal
5. menthone
6. citronellal
7. methyl nonanoate
8. 2,3-butanediol
9. linalool
10. linalyl acetate
11. methyl decanoate
12. menthol
13. α -terpineol
14. γ -terpineol
15. methyl undecanoate

GC_FF00433

50m, 0.32mm ID, 0.33 μ m, traditional Carbowax® 20M or Rt-CW20M™ F&F (cat.# 12539)
On-column injection of 5ng to 150ng each compound in methylene chloride, split 10:1

Carrier gas: hydrogen, 40cm/sec.
Inj. & det. temp.: 220°C
Oven temp.: 110°C isothermal

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com