


Alcoholic Beverages Acids and Esters

Caproic Acid and Ethyl Laurate

Stabilwax®-DA

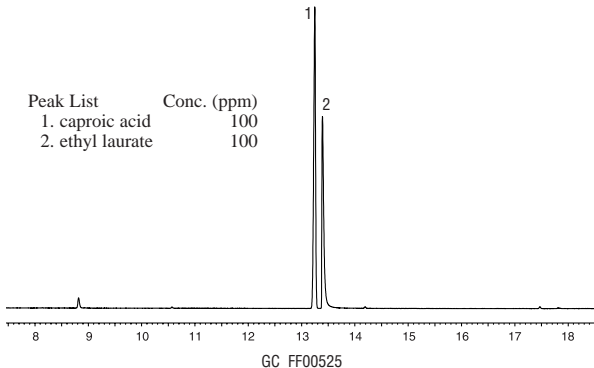

Stabilwax®-DA 30m, 0.18mm ID, 0.18 μ m (cat.# 550752)

Inj.: 1 μ L splitless (hold 0.5 min.) at conc. shown in peak list, in ethyl acetate, 4mm ID splitless liner w/wool (cat.# 20814-202.1)

Inj. temp.: 240°C

Carrier gas: hydrogen

Make-up gas: nitrogen

Linear velocity: 28psi @ 240°C

Oven temp.: 80°C to 230°C at 5°C/min.

Det.: FID