


5% Fragrance Materials Association Mix Rtx[®]-1

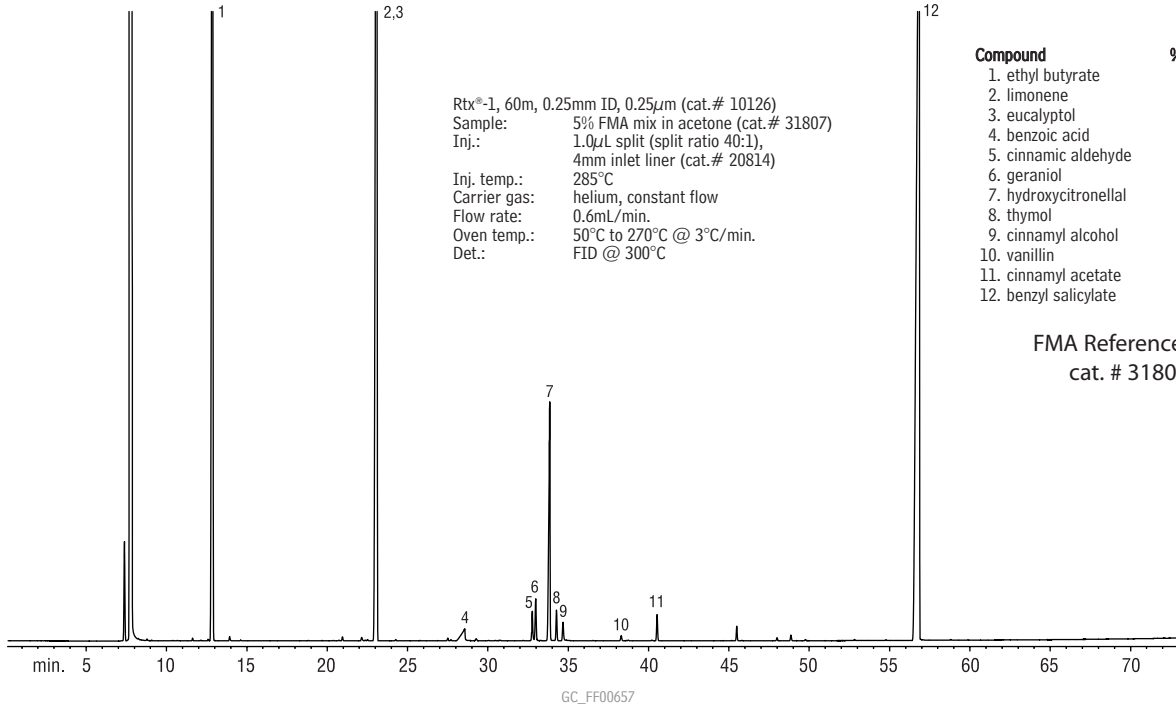

Rtx[®]-1, 60m, 0.25mm ID, 0.25 μ m (cat.# 10126)
Sample: 5% FMA mix in acetone (cat.# 31807)
Inj.: 1.0 μ L split (split ratio 40:1),
4mm inlet liner (cat.# 20814)
Inj. temp.: 285 $^{\circ}$ C
Carrier gas: helium, constant flow
Flow rate: 0.6mL/min.
Oven temp.: 50 $^{\circ}$ C to 270 $^{\circ}$ C @ 3 $^{\circ}$ C/min.
Det.: FID @ 300 $^{\circ}$ C

Compound	% in the FMA mix
1. ethyl butyrate	35.7
2. limonene	20.0
3. eucalyptol	0.5
4. benzoic acid	1.0
5. cinnamic aldehyde	0.5
6. geraniol	0.6
7. hydroxycitronellal	5.0
8. thymol	0.3
9. cinnamyl alcohol	0.3
10. vanillin	0.1
11. cinnamyl acetate	0.3
12. benzyl salicylate	35.7

FMA Reference Mix
cat. # 31807