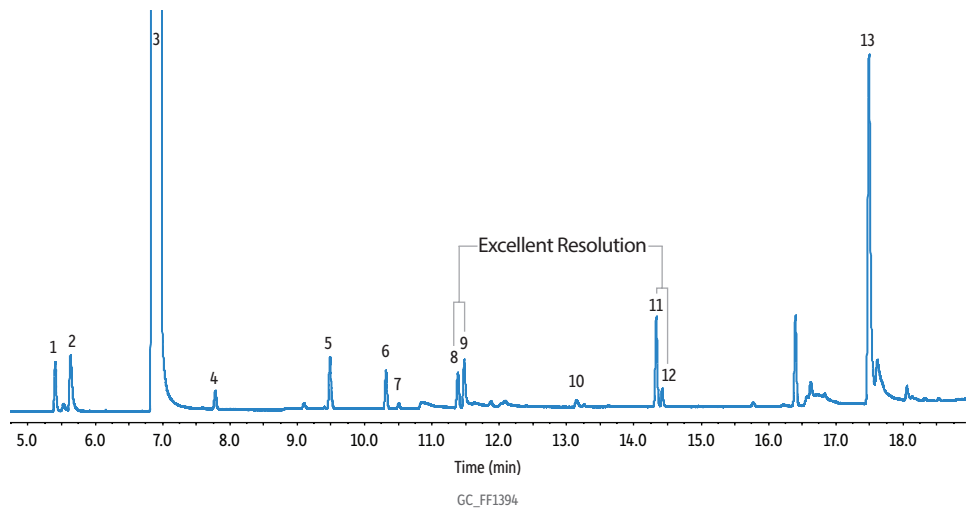




Peaks	tr (min)
1. Acetaldehyde	5.404
2. Methanol	5.632
3. Ethanol	6.956
4. Ethyl formate	7.779
5. 1-Propanol	9.485
6. Ethyl acetate	10.318
7. sec-Butanol	10.507
8. Isobutanol	11.384
9. Acetic acid	11.479
10. Acetal	13.143
11. 3-Methyl-1-butanol	14.333
12. 2-Methyl-1-butanol	14.419
13. Furfural	17.491

Column	Rtx-1301, 60 m, 0.25 mm ID, 1.0 µm (cat.# 16056)
Standard/Sample	Marsala (neat)
Diluent:	n/a
Injection	
Inj. Vol.:	1 µL split (split ratio 100:1)
Liner:	Topaz precision inlet liner, 4.0 mm x 6.3 x 78.5, with wool (cat.# 23305)
Inj. Temp.:	250 °C
Split Vent Flow Rate:	144.4 mL/min
Oven	
Oven Temp.:	35 °C (hold 3 min) to 220 °C at 8 °C/min (hold 5 min)
Carrier Gas	He, constant flow
Flow Rate:	1.4 mL/min
Linear Velocity:	25.457 cm/sec @ 35 °C
Detector	FID @ 250 °C
Make-up Gas Flow Rate:	25 mL/min
Make-up Gas Type:	N ₂
Hydrogen flow:	35 mL/min
Air flow:	400 mL/min
Data Rate:	50 Hz
Instrument	Agilent 7890A GC