



Unleaded Gasoline Rtx®-1

Unweathered

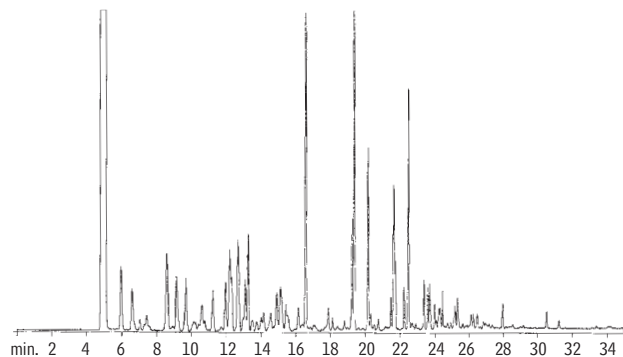

99% Weathered

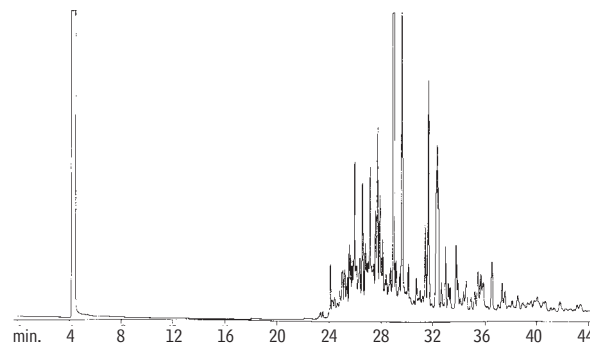

GC_MS00270

30m, 0.53mm ID, 1.50 μ m Rtx®-1 (cat.# 10170)
Oven temp.: 40°C (hold 3 min.) to 75°C @
15°C/min. to 275°C @ 20°C/min. (hold 5 min.)
Inj./det. temp.: 250°C/285°C
Carrier gas: hydrogen
Linear velocity: 50cm/sec. set @ 40°C
FID sensitivity: 4.10 x 10⁻⁹ AFS
Split ratio: 30:1

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com