


## Permanent Gases Rt-MSieve™ 5A PLOT

1. hydrogen, 29.1 ppm
2. argon, 53.4 ppm
3. oxygen, 31.3 ppm

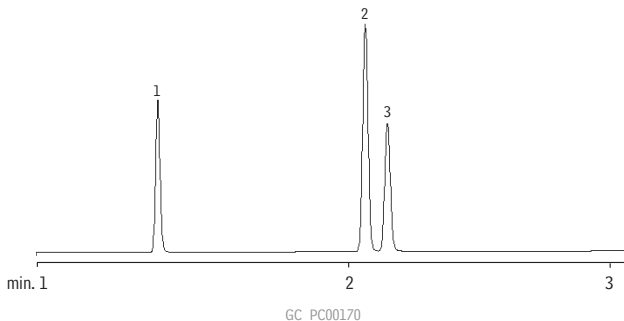

30m, 0.53mm ID, 50 $\mu$ m Rt-MSieve™ 5A PLOT (cat.# 19723)

Sample: 0.5mL (sample loop injection)

Column temp.: 27°C  
Carrier gas: helium  
Linear velocity: 34cm/sec.  
Det.: Valco HID

Chromatogram courtesy of Larry McElmurry, Mobile Analytical Labs.

Restek Corporation 110 Benner Circle Bellefonte, PA 16823  
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • [www.restek.com](http://www.restek.com)