


Hydrocarbon Gases

HayeSep® S

(micropacked)

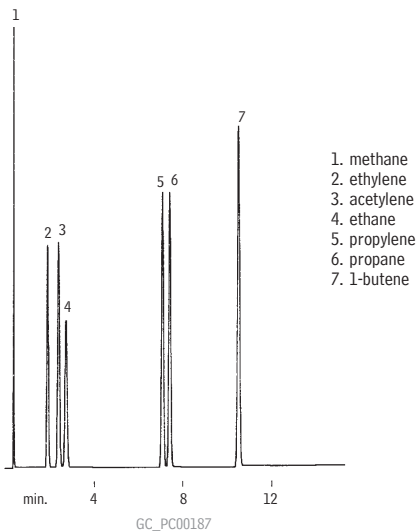

2m, 0.75mm ID HayeSep S micropacked column (cat.# 19011)
500 μ L split injection of a light hydrocarbon gas mixture

Oven temp.: 40°C (hold 3 min.) to 150°C @
15°C/min. (hold 5 min.)
Inj. & det. temp.: 220°C
Carrier gas: helium
Flow: 20mL/min. set @ 40°C
FID sensitivity: 32 x 10⁻¹⁴ AFS
Split ratio: 10:1

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com