


Aromatics (Benzene/Toluene/Xylene) Rtx[®]-200

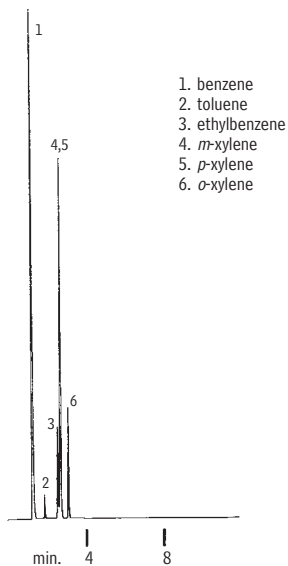

GC_PC00189

30m, 0.53mm ID, 0.50 μ m Rtx[®]-200 (cat.# 15040)
0.1 μ L direct injection of BTEX Standard (cat.# 30051)

Oven temp.: 60°C
Inj. & det. temp.: 250°C
Carrier gas: hydrogen
Linear velocity: 40cm/sec. (flow rate: 5.2cc/min.)
FID sensitivity: 4 x 10⁻¹¹ AFS