


Aromatics (Benzene/Toluene/Xylene) Rtx[®]-Wax

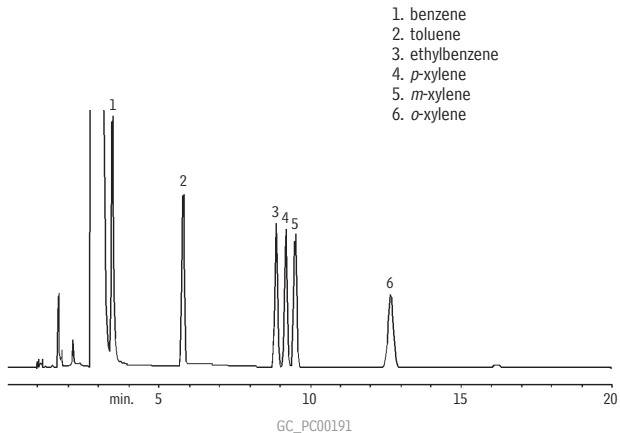

30m, 0.53mm ID, 1.0 μ m Rtx[®]-Wax (cat.# 12455)
0.1 μ L direct injection of BTEX Standard (cat.# 30051)
Oven temp.: 50°C
Inj. & det. temp.: 250°C
Carrier gas: hydrogen
Linear velocity: 40cm/sec. (flow rate: 5.2cc/min.)
FID sensitivity: 16 x 10⁻¹¹ AFS