


Propane Fuel Rtx®-1

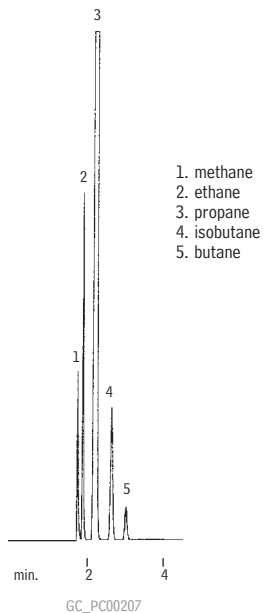

60m, 0.53mm ID, 5.0 μ m Rtx®-1 (cat. # 10183)

1.0 μ L propane fuel

Direct injection using a Uniliner® inlet liner

Oven temp.: 40°C

Inj. & det. temp.: 200°C

Carrier gas: hydrogen

Linear velocity: 80cm/sec. (flow rate: 10cc/min.)

FID sensitivity: 4 x 10⁻¹¹ AFS

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com