


# Lubricant (Synthetic) Stabilwax®

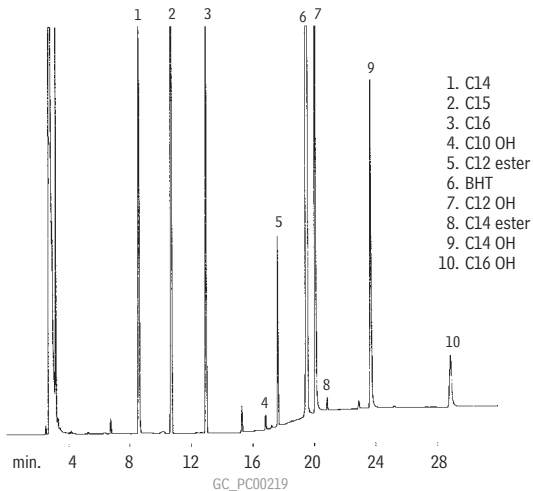

30m, 0.53mm ID, 1.0 $\mu$ m Stabilwax® (cat.# 10655)  
1.0 $\mu$ L split injection of a synthetic lubricant  
Oven temp.: 100°C to 175°C @ 5°C/min., to 225°C  
@ 15°C/min. (hold 10 min.)  
Inj. / det. temp.: 250°C / 275°C  
Carrier gas: helium  
Linear velocity: 20cm/sec. (flow rate: 2.6cc/min.)  
FID sensitivity: 8 x 10<sup>-11</sup> AFS  
Split vent: 40:1