


Simulated Distillation (High-Temperature Calibration) MXT[®]-1 Sim Dist

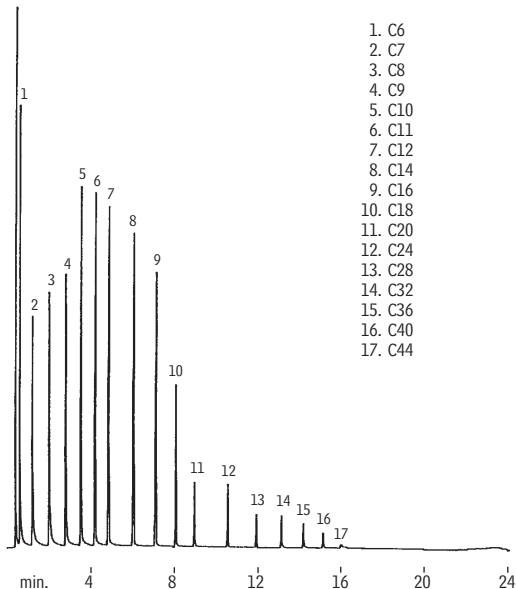

GC_PC00224

6m, 0.53mm ID, 0.15 μ m MXT[®]-1 Sim Dist (cat.# 70101)
Wet needle on-column injection of ASTM D2887 standard.
Oven temp.: -12°C to 430°C @ 20°C/min.
Inj. & det. temp.: -17°C to 433°C/430°C
Carrier gas: hydrogen
Linear velocity: 40cm/sec.
FID sensitivity: 128 x 10⁻¹¹ AFS