


# Crude Oil, Arabian Light Rtx<sup>®</sup>-1

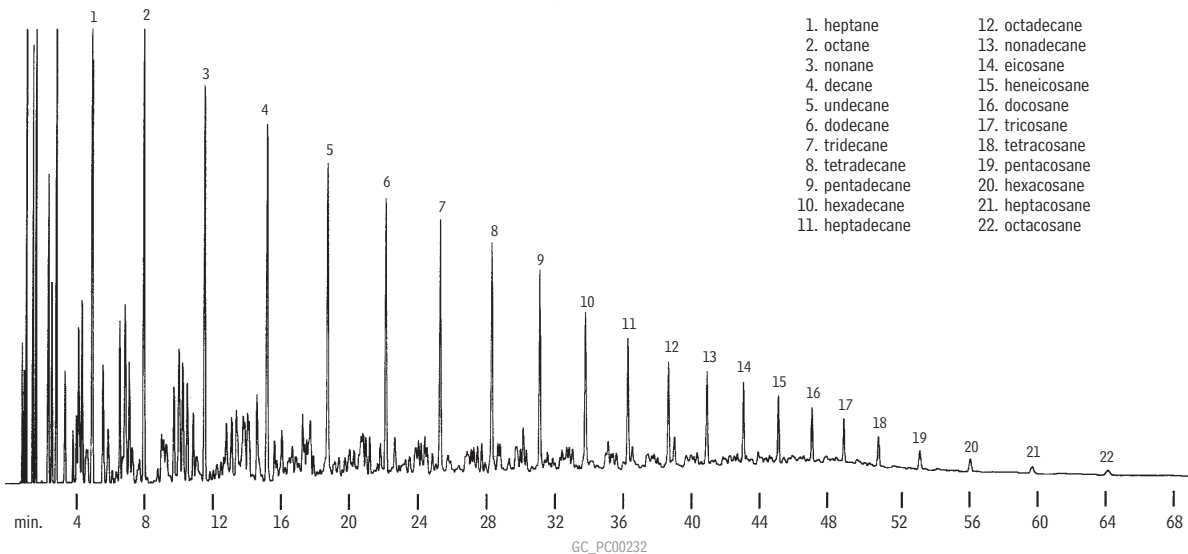

30m, 0.53mm ID, 3.0 $\mu$ m Rtx<sup>®</sup>-1 (cat.# 10185)  
0.6 $\mu$ L split injection of an Arabian light crude oil.  
Oven temp.: 40°C to 285°C @ 5°C/min. (hold 15 min.)  
Inj. & det. temp.: 285°C  
Carrier gas: hydrogen  
Linear velocity: 40cm/sec. (flow rate: 5.2cc/min.)  
FID sensitivity: 128 x 10<sup>-14</sup> AFS  
Split ratio: 50:1