


OPN on Res-Sil™ C packing demonstrates unique selectivity for *cis*-2-butene and 1,3-butadiene.

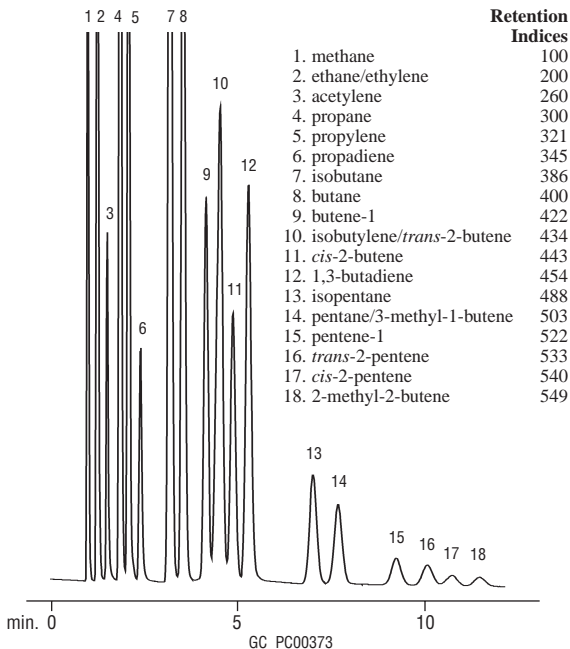

OPN on Res-Sil™ C Packing, 80/100 mesh, 12' x 2mm ID x 1/8" OD in SilcoSmooth™ tubing (cat.# 80437). 20µL on-column injection of refinery gas.

Concentration: 0.1-6 absolute mole %
 Oven temp.: 50°C
 Inj. & det. temp.: 200°C
 Carrier gas: He
 Flow rate: 30mL/min

Standard courtesy of AC Analytical Controls, Bensalem, PA.

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restekcorp.com