


A refinery gas calibration standard on a Restek refinery gas packed column set.

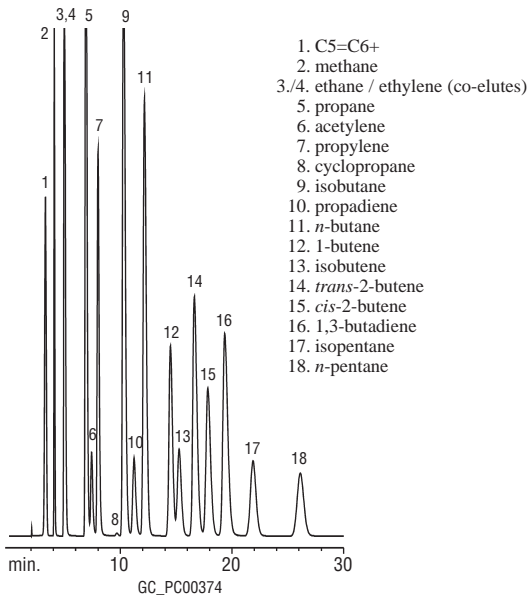

Refinery Gas Column Set. This is a 3 column set identified as 2abc Refinery Gas Set

Oven temp.: 60°C
Inj. temp.: 150°C
Det. temp.: 150°C FID
Flow rate: 30mL/min., helium
Sample: refinery gas C1-C5
Sample size: 1cc

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restekcorp.com