


Natural Gas ShinCarbon ST (Micropacked)

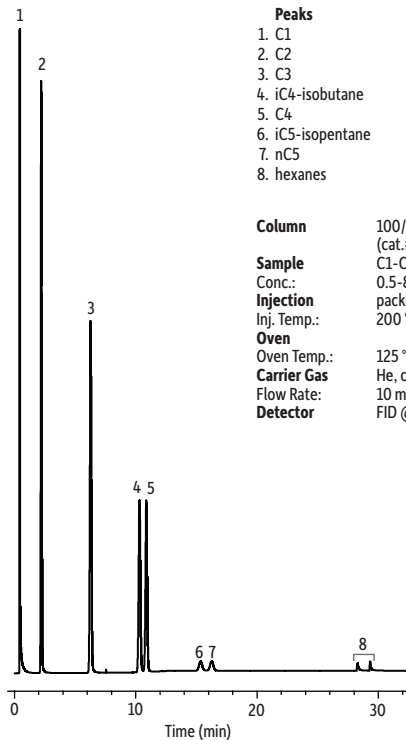

- Peaks**
1. C1
 2. C2
 3. C3
 4. iC4-isobutane
 5. C4
 6. iC5-isopentane
 7. nC5
 8. hexanes

Column 100/120 mesh on Shincarbon ST, SilcoSmooth tubing, 1 m, 1/16 OD, 1.0 mm ID (cat.# 19809)
Sample C1-C6 hydrocarbons
Conc.: 0.5-8% varied
Injection packed on-column
Inj. Temp.: 200 °C
Oven
Oven Temp.: 125 °C to 250 °C at 16 °C/min (hold 25.2 min)
Carrier Gas He, constant flow
Flow Rate: 10 mL/min
Detector FID @ 250 °C

GC_PC00672



© 2020 Restek Corporation. All rights reserved.
www.restek.com

Restek patents and trademarks are the property of Restek Corporation. (See www.restek.com/Patents-Trademarks for full list.) Other trademarks in Restek literature or on its website are the property of their respective owners. Restek registered trademarks are registered in the U.S. and may also be registered in other countries.