


Ketones

Rt®-QS-BOND

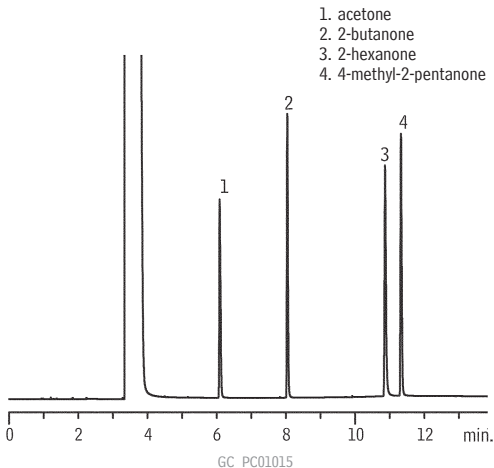

Column: Rt®-QS-BOND, 30m, 0.53mm ID (cat. # 19738)
Sample: 1,000µg/mL VOA Calibration Mix #1 (ketones)
(cat. # 30006) in methanol:water, 90:10 (v/v)
Inj.: 1µL split (split ratio 50:1), 4mm Siltek® single
goose-neck liner (cat. # 20798-214.1)
Inj. temp.: 240°C
Carrier gas: helium, constant flow
Flow rate: 5.7mL/min. @ 40°C
Oven temp.: 60°C to 225°C @ 15°C/min. (hold 10 min.)
Det.: FID @ 240°C