


Halogenated Hydrocarbons on Rt®-Alumina BOND/CFC

Peaks

- 1,1-Dichloroethylene
- Methylene chloride
- trans*-1,2-Dichloroethylene
- Tetrachloromethane
- Unknown
- cis*-1,2-Dichloroethylene
- Trichloroethylene
- 1,1,1-Trichloroethane
- Chloroform

Column Rt®-Alumina BOND/CFC, 30 m,
0.53 mm ID, 10 μ m (cat.# 19763)
Injection Split Vent
Flow Rate: 170 mL/min.
Oven
Oven Temp: 135 °C
Carrier Gas H₂, constant flow
Flow Rate: 8 mL/min. @ 135 °C
Detector FID @ 220 °C

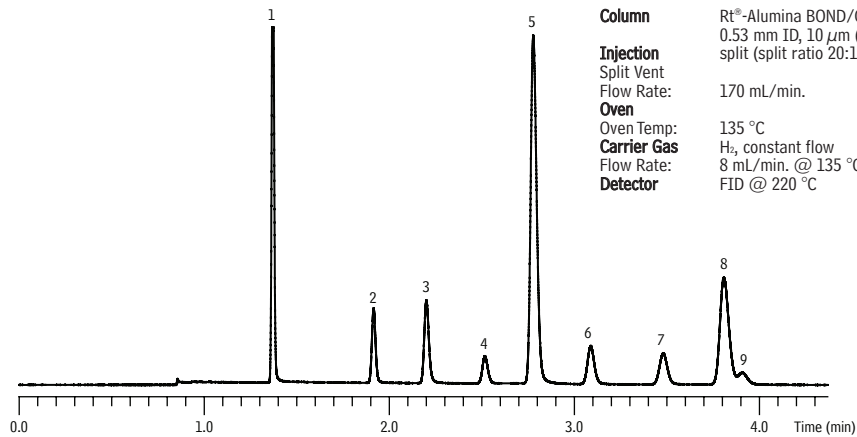

GC_PC1193