


Sulfur Compounds on Rt[®]-XLSulfur (micropacked)

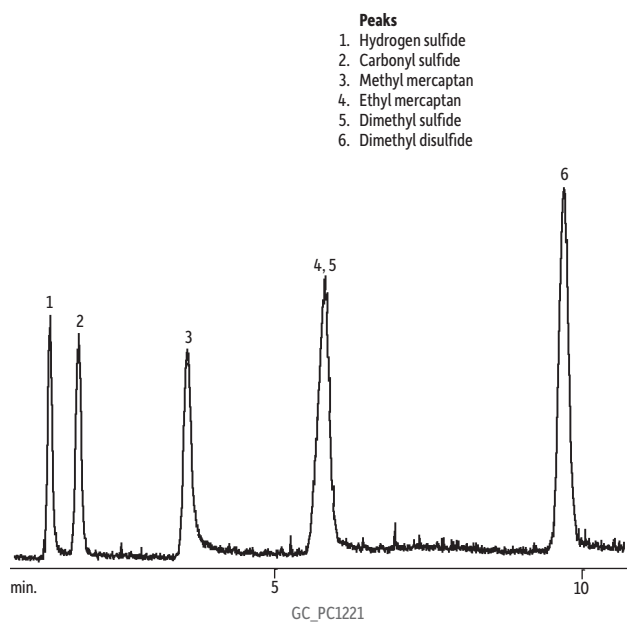

| | |
|--------------------|--|
| Column | Rt [®] -XLSulfur, 1 m, 0.75 mm ID (cat.# 19806) |
| Sample | |
| Conc.: | 1 mL of 50 ppbv each sulfur compound |
| Injection | sample valve |
| Oven | |
| Oven Temp: | 60 °C to 230 °C at 15 °C/min. |
| Carrier Gas | He, constant flow |
| Flow Rate: | 9 mL/min. |
| Detector | SCD |

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com