


Alkaloids on Diatomaceous Earth WHP Packed Column

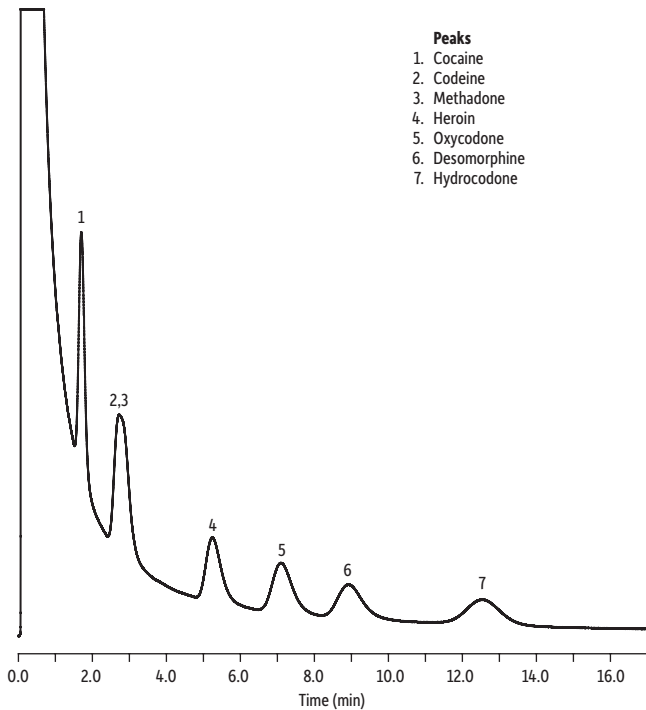

GC_PC1296

- Peaks**
1. Cocaine
 2. Codeine
 3. Methadone
 4. Heroin
 5. Oxycodone
 6. Desomorphine
 7. Hydrocodone

Column	3% OV-17, 0.1% Versamid 900, 100/120 mesh on diatomaceous earth WHP, SilcoSmooth tubing, 1 m, 1/8 in OD, 2 mm ID
Sample	Methadone (cat.# 34005) Cocaine (cat.# 34015) Codeine (cat.# 34000) Morphine (cat.# 34006) 140 µg/mL
Conc.:	
Injection	
Inj. Vol.:	1 µL on-column
Inj. Temp.:	250 °C
Oven	
Oven Temp.:	230 °C (hold 15 min)
Carrier Gas	He, constant flow
Flow Rate:	40 mL/min
Detector	FID @ 275 °C
Instrument	Agilent 7890B GC