


Opiates (TMS Derivatives) Rtx®-5

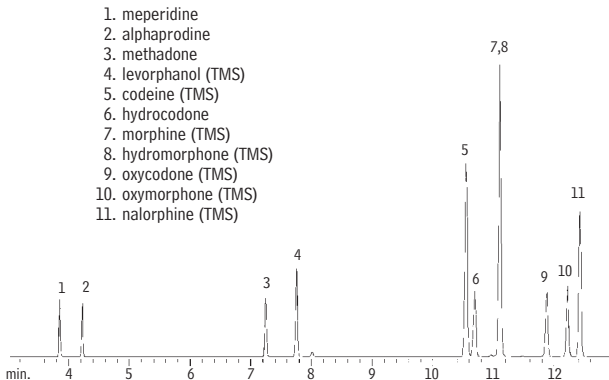

GC_PH00233

30m, 0.25mm ID, 0.25 μ m Rtx®-5 (cat.# 10223)
2.0 μ L split injection of opiates
Conc: 2000ng/ μ L
Oven temp.: 200°C to 325°C @ 7°C/min.
Inj. / det. temp.: 250°C / 300°C
Det. type: MS
Ionization: EI
Carrier gas: helium
Mode: full scan
Linear velocity: 30cm/sec. set @ 200°C
Split ratio: 50:1