


Phencyclidine Rtx®-5

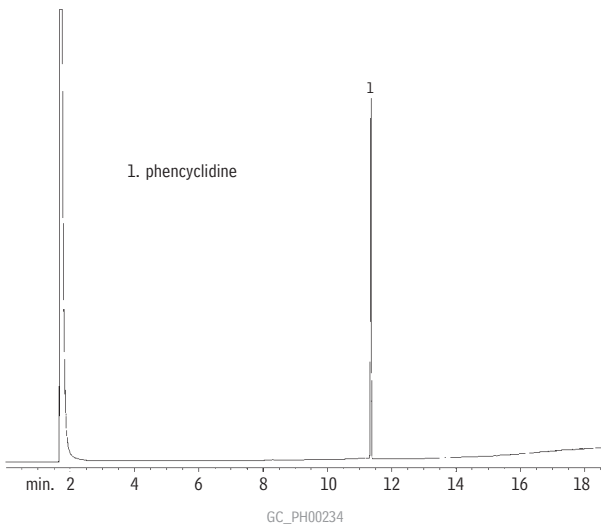

30m, 0.25mm ID, 0.25 μ m Rtx®-5 (cat.# 10223)
1.0 μ L split injection of phencyclidine
Conc.: 1000ng/ μ L
Oven temp.: 50°C (hold 1 min.) to 250°C @ 25°C/min.,
then to 325°C @ 10°C/min. (hold 2 min.)
Inj. / det. temp.: 250°C / 325°C
Carrier gas: helium
Linear velocity: 30cm/sec. set @ 50°C
FID sensitivity: 2.56 x 10¹⁰ AFS
Split ratio: 30:1