


# Fenfluramine (TFAA Derivative) Rt- $\beta$ DEXcst™

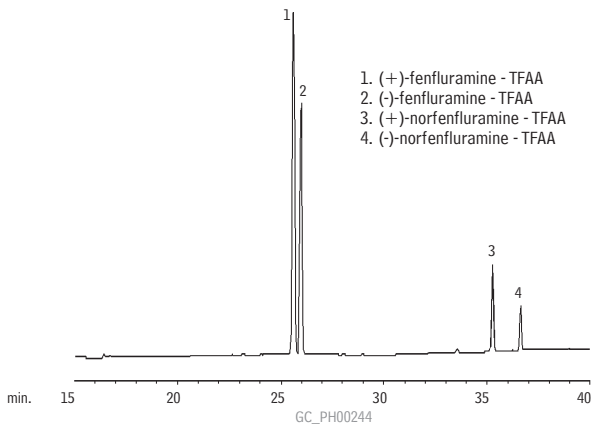

30m, 0.32mm ID, 0.25 $\mu$ m  
Rt- $\beta$ DEXcst™ (cat.# 13102)  
On-column injection, 25ng/enantiomer

Oven temp.: 90°C (hold 1 min.) to 200°C  
@ 2°C/min. (hold 3 min.)  
Carrier gas: hydrogen  
Linear velocity: 80cm/sec. set @ 60°C  
Det.: FID set @ 220°C