


## Cocaine & Metabolites (TMS Derivatives) Rtx®-5

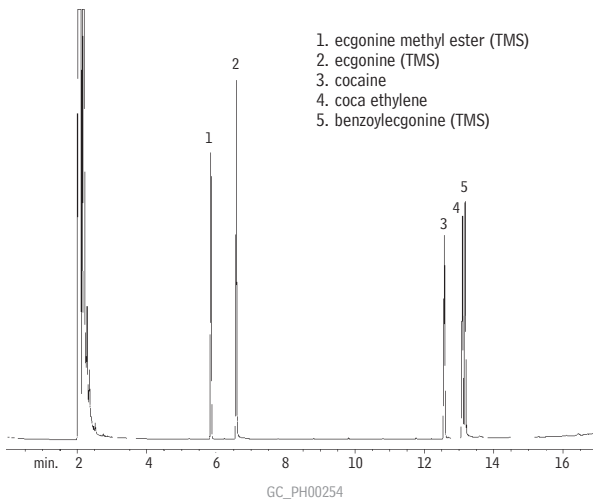

30m, 0.25mm ID, 0.25 $\mu$ m Rtx®-5 (cat.# 10223)  
1.0 $\mu$ L split injection of cocaine and cocaine metabolites

Oven temp.: 150°C to 320°C @ 10°C/min.  
Inj. / det. temp.: 250°C / 300°C  
Carrier gas: helium  
Linear velocity: 30cm/sec. set @ 50°C  
FID sensitivity: 2.56 x 10<sup>-10</sup> AFS  
Split ratio: 30:1

Restek Corporation 110 Benner Circle Bellefonte, PA 16823  
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com