


Anesthetics

Rtx[®]-BAC1 & Rtx[®]-BAC2


(headspace analysis)

1. isoflurane
2. enflurane
3. halothane
4. methoxyflurane

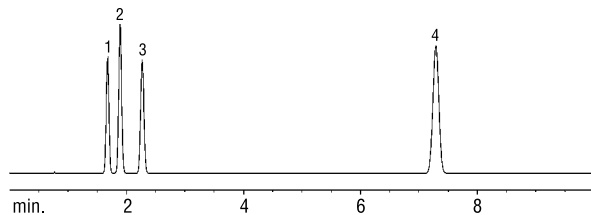

GC_PH00501

30m, 0.53mm ID, 3.0µm Rtx[®]-BAC1 (cat.# 18001)

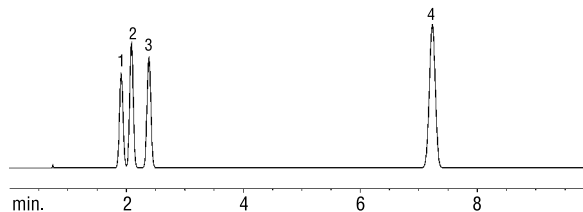

GC_PH00502

30m, 0.53mm ID, 2.0µm Rtx[®]-BAC2 (cat.# 18000)

1.0mL headspace sample
Oven temp.: 40°C (hold 5 min.) to 240°C @ 5°C/min.
Inj. & det. temp.: 240°C
Carrier gas: helium
Linear velocity: 65cm/sec.