


Antibiotic Amoxicillin on Ultra IBD

Peak

1. Amoxicillin

Column Ultra IBD (cat.# 9175565)
Dimensions: 150 mm x 4.6 mm ID
Particle Size: 5 μm
Pore Size: 100 \AA
Temp.: 30 $^{\circ}\text{C}$
Sample
Diluent: acetonitrile:water (10:90, v/v)
Conc.: 1.5 mg/mL
Inj. Vol.: 5 μL
Mobile Phase 10 mM ammonium formate, pH 2.5:acetonitrile (95:5)
Flow: 1.2 mL/min
Detector UV/Vis @ 270 nm

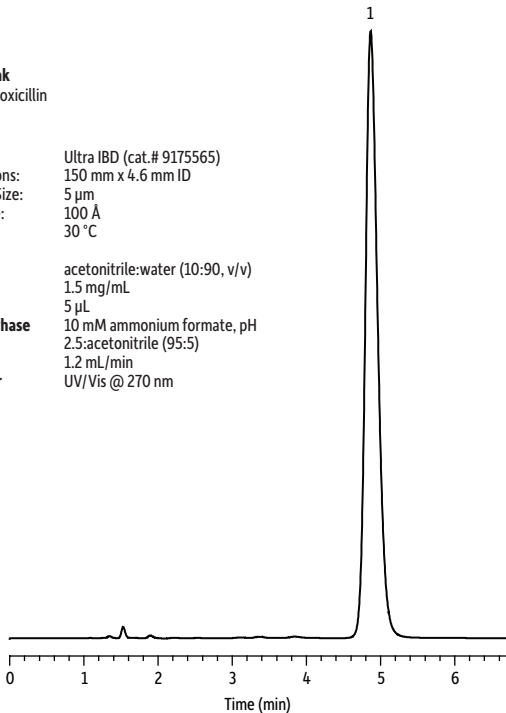

LC_0095