


## Antibiotic Cephaloridine on Ultra IBD (Normal Phase)

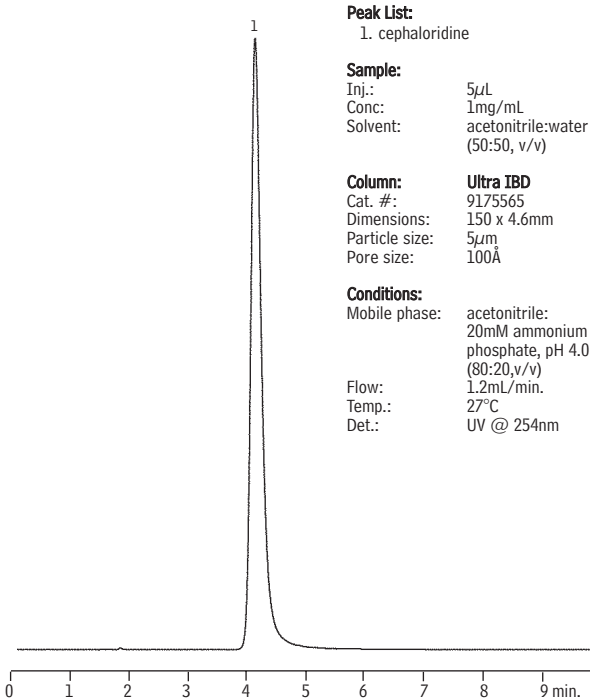

LC\_0101

### Peak List:

1. cephaloridine

### Sample:

Inj.: 5 $\mu$ L  
Conc: 1mg/mL  
Solvent: acetonitrile:water  
(50:50, v/v)

### Column: Ultra IBD

Cat. #: 9175565  
Dimensions: 150 x 4.6mm  
Particle size: 5 $\mu$ m  
Pore size: 100Å

### Conditions:

Mobile phase: acetonitrile:  
20mM ammonium  
phosphate, pH 4.0  
(80:20,v/v)  
Flow: 1.2mL/min.  
Temp.: 27°C  
Det.: UV @ 254nm